

**>> OUR CONVENIENT  
AUTOCLAVE, RIGHT ON  
THE LABORATORY  
BENCH. <<**



**18 LITER VOLUME**

**STERILISATION**

**MATERIALS TESTING**

**MEDIA PREPARATION**


# Applications

## THE WIDE RANGE OF USING CERTOCLAV AUTOCLAVES



### MICROBIOLOGY / BIO-CHEMISTRY

- Autoclave culture mediums
- Refuse sterilisation
- Efficiently processes small amounts – just in time

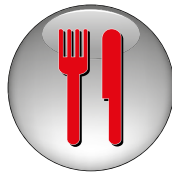

### FOOD INDUSTRY

- The complete autoclave solution for laboratories in the food processing industry
- Sterilisation of utensils and refuse
- Culture medium preparation



### MATERIAL TESTING

- Testing for vapour resistance and material quality
- Long running cycles, ICP-TM 650, HAST test, pressure vessel test, hydrolysis, GLB-AL 631
- Special fabrications available upon request
- CertoClav MultiControl with a freely programmable electronic control



### SCHOOLS AND UNIVERSITIES

- Robust autoclave, ideal for daily training operation
- Special deals for schools – please inquire



### PHARMACY

Sterilisation of instruments and packaging (i.e. eyedrop bottles and liquids)

# Overview

## THE MAIN FEATURES OF OUR PRODUCTS




		Max. Temperature	Electronic Control	Programmable	Touch Screen + WLAN	USB Port	Printer Port	Volume	Chamber material	
ESSENTIAL	121	✓					✓	12/18	Alu	<ul style="list-style-type: none"> <li>Affordable solution</li> <li>Printer for documentation</li> <li>Automatic</li> </ul>
CLASSIC	140							18	SS	<ul style="list-style-type: none"> <li>Very low maintenance (no electronics)</li> <li>Time tested for 50 years</li> </ul>
MULTICONTROL 2	140	✓	✓		✓	✓		18	SS	<ul style="list-style-type: none"> <li>Programmable autoclave</li> <li>Multi-step time/temperature</li> </ul>
MULTICONTROL 2 TOUCH	140	✓	✓	✓	✓	✓		18	SS	<ul style="list-style-type: none"> <li>Same as MultiControl 2</li> <li>Touch Screen Controller</li> <li>Cloud access</li> </ul>
CONNECT	140	✓	✓	✓		✓*		18	SS	<ul style="list-style-type: none"> <li>Fully automatic (electric valves)</li> <li>Media sensor optional</li> </ul>
CULTURE MEDIA	140	✓	✓	✓		✓*		18**	SS	<ul style="list-style-type: none"> <li>Media preparator</li> <li>Internal stirrer</li> <li>Port for pumping, sterile Medien</li> </ul>

\* via Tablet Controller    \*\* 10 liter Media Container




# CertoControl

## THE USER INTERFACE WITH ON-SCREEN TRAINING.

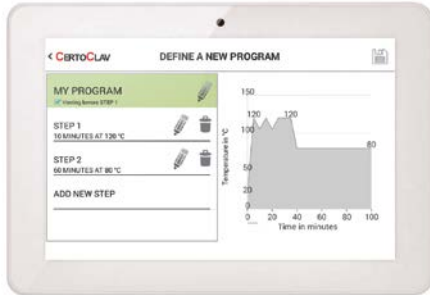

### Define custom programs

- With a simple interface you can program multi-step processes of temperature and time. E.g. 20 minutes at 110° C, then 3 hours at 120° C and then 1 hours at 80° C.
- Long term cycles possible. Programming can be done in minutes, hours and days.
- Program is graphically displayed.
- Programs can be named.
- One-click solution to program a fractioned venting process.

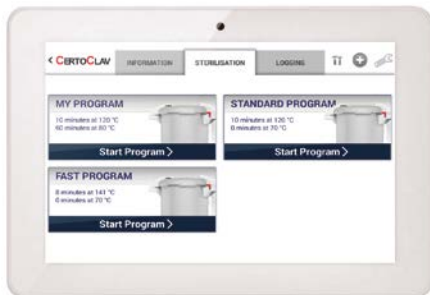


### Choose program

- All programs are displayed as icons with name and key data.
- Programs can be added and deleted.
- Full list of steps of each program can be displayed.



### Watch step-by-step video instructions

- All usage steps are explained in videos, e.g. when used as an autoclave the steps to load the unit are explained.
- No more need to train staff – the unit explains everything itself.
- Videos can be skipped once staff is proficient in using the equipment.

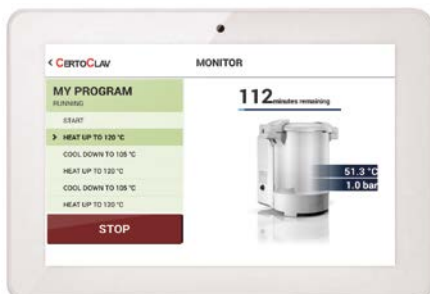

### Show status of running programs

- Current step of program is displayed.
- Current data displayed.
- Time passed and remaining is displayed.



### Scroll through protocol history

- Protocols are stored on the device.
- Target temperature and actual temperature are graphically displayed and accessible in numerical format.
- Error messages are displayed in detailed with explanation.

# CertoCloud

EASY DATA ACCESS AND STORAGE.  
ONLINE DEVICE MANAGEMENT.



## Secure Connection

- SSL encryption
- Secure server
- 100 % data privacy




## Multi Device Management

- Combine all devices in one account
- Create programs
- Disseminate programs to all devices



## Long Term Data Storage

- All data stored for 10 years to fulfil documentation requirements
- Secure backup

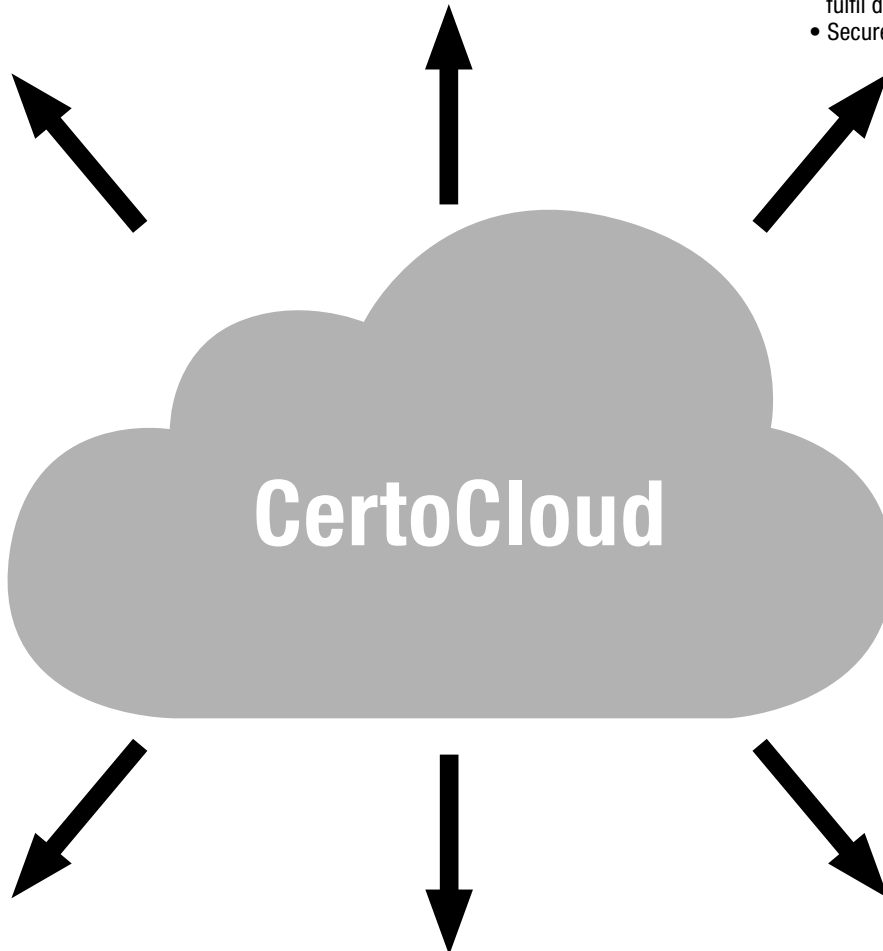

CertoCloud



## Customization

- Custom reports
- Direct printing
- Freely definable data formats



## Email Notifications

- Alarms
- Cycle complete
- Process data



## User Management

- Individual accounts
- Accounts linked to devices



# CertoClav Essential

## THE ECONOMIC BASIC AUTOCLAVE

LOW COST  
SINGLE BUTTON  
CONTROL  
FRACTIONED  
VENTING



Art.-No. 8510173 (12 litres) / 8510183 (18 litres)

### FEATURES

- Autoclave with one button controller. Just push „Start“ and the sterilisation process automatically runs.
- Sterilisation process consists of:
  - 3 x reaching 115° C and cooling for fractioned venting
  - Reaching 121° C and holding for 15 Minutes
  - Cool down to 80° C
- Printer port for printing protocols.

### INCLUDED

- Control unit with single button control
- Printer port

### ACCESSORY

#### Wire basket

Art.-No. 8583010  
Diam. 23, stainless steel  
For easy filling with  
sterilisation goods.



#### Printer

Art.-No. 1300187  
For easy print out of  
your documentation.



# CertoClav Classic

**AUTOCLAVING UP TO 140°C  
WITH TIME TESTED TECHNOLOGY**

**THERMOSTAT  
PRECISION VALVE  
TIMER**



Art.-No. 8510174 (for valve with 125 / 140° C), 8510184 (for valve with 115 / 121° C). Both 18 litres.

## FEATURES

- One of the lowest cost solutions for reliable sterilisation and materials testing up to 140° C / 2,7 bar.
- Ultra-low maintenance autoclave with time tested technology.
- Precision valves guarantee that the required pressure and temperature for sterilisation are reached.
- Two levels selectable:  
125° C (1.4 bar / 20 psi) or 140° C (2.7 bar / 39 psi),  
alternative valve for 115° C / 121° C available

## INCLUDED

- Precision valve for 125 / 140° C (intermediate values possible)
- Thermostat for heat regulation
- Timer switch (5 - 60 minutes, permanent on)

## ACCESSORY

### Wire basket

Art.-No. 8583010  
Diam. 23, stainless steel

For easy filling with  
sterilisation goods.



# CertoClav MultiControl 2

## THE FREELY PROGRAMMABLE AUTOCLAVE

PRINTER PORT  
MULTI-LANGUAGE  
PC PORT



Art.-No. 8510175 (18 litres)

### FEATURES

- Set programs for all sterilisation tasks
- Logging of processes with PC Software (basic version included) or via printer.
- Simple program programmable on the controller (1 time / temperature step)
- Complex programs (multiple time / temperature steps) can be entered with PC Software (Professional version required).
- Automatic venting process for reliable sterilisation.
- Programmable from room temperature to 140° C
- Suitable for materials testing (ICP-TM610, HAST, GLB-AL 631)

### INCLUDED

- Controller with with printer port and USB port
- Large LCD screen with temperature display

### ACCESSORY

#### Wire basket

Art.-No. 8583010  
Diam. 23, stainless steel  
For easy filling with  
sterilisation goods.



#### Printer

Art.-No. 1300187  
For easy print out of  
your documentation.



#### PC Software Professional Version

Art.-No. 8200210 (Art.-No. 8200211  
when purchased in bundle with autoclave)  
For programming of multi step  
programs.



# CertoClav MultiControl 2 Touch

## THE AUTOCLAVE WITH TOUCH SCREEN AND WIFI CONNECTION

PROGRAMMABLE

TOUCH SCREEN  
CONTROLLER

WLAN/CLOUD  
CONNECTION



Art.-No. 8510298 (18 litres)

### FEATURES

- Automatic long term documentation in the cloud. Digitally signed protocols according to 21 CFR Part/GMP guidelines.
- Freely programmable autoclave (by time and temperature, 6 program slots, each 14 steps).
- All features of CertoControl and CertoCloud.
- Printer port for printing of local.

### INCLUDED

- 7 inch HD touch screen controller with WiFi and Bluetooth.
- Software with key features:
  - User Management
  - On-Screen video training and instructions for use
  - Live update of software
  - Video connection to support
  - One-touch ordering of spareparts
  - Programming of autoclave
  - Multiple languages
- Additional stand-alone power electronics that can be operated without touch screen controller.

### ACCESSORY

#### Wire basket

Art.-No. 8583010  
Diam. 23, stainless steel  
For easy filling with  
sterilisation goods.



#### Printer

Art.-No. 1300187  
For easy print out of  
your documentation.



#### Tablet controller holder

Art.-No. 1300463  
Wall and table mountable  
holding arm.



# CertoClav Culture Media

## THE AUTOCLAVE FOR EASY MEDIA PREPARATION

MEDIA  
PREPARATOR

CFR 21 PART 11  
COMPLIANT CLOUD  
DOCUMENTATION

PROCESS PORTS



Art.-No. 8510177 (18 litres)

### FEATURES

- Media preparator for up to 10 litres of culture media. Sterilisation up to 140° C.
- Ports for pumping out liquid and inserting additives.
- In waterbath mode, temperatures under 100° C can be operated with open lid.

### INCLUDED

- 10 litre culture media stainless steel bucket.
- Media sensor for process guidance based on media temperature.
- Magnetic Stirrer for homogeneity.
- Process ports for additives and for pumping out sterile media.
- Touch screen controller with cloud connection for GMP / 21 CFR Part 11 compliant documentation (3 month cloud usage included).

### ACCESSORY

#### Wire basket

Art.-No. 8583010  
Diam. 23, stainless steel  
For easy filling with  
sterilisation goods.



#### Tablet controller holder

Art.-No. 1300463  
Wall and table mountable  
holding arm.

# TECHNICAL DETAILS

	ESSENTIAL	CLASSIC	MULTICONTROL 2
Key benefit	Low cost professional autoclave	Low cost autoclave for sterilisation up to 140° C with low maintenace	Electronic autoclave with PC/ printer documentation. Freely programmable (with Pro Software)
Core use case	Automatic sterilisation of solids and waste	Sterilisation with low requirement concerning documentation	Sterilisation up to 140° C with documentation on printer or PC, Materials testing
Temperature control	Electronic	Mechanical thermostat	Electronic
Pressure control	Via temperature	Mechanical precision valve	Via temperature
Interface	Single button and 5 LEDs	Temperature selection knob	LCD display
Steam release	Automatic Steam release cock	Manual steam release cock	Automatic Steam release cock
Other key components	Printer port	Mechanical timer	
Optional equipment	Wire basket, printer	Wire basket	Wire basket, printer, pro licence of PC software (for multi step programs)
Temperature range	121° C	125° C / 140° C selectable. Altern. model: 115° C / 121° C	Room temp. to 140°C, water bath mode with open lid
Max working pressure	1,2 bar	2,7 bar	
Chamber material	Aluminium	Stainless steel 304L	
Volume	12 or 18	18	
Diameter	24	26	
Rated voltage, Output	230V AC, 50-60 Hz, 1900W	230 V AC, 50-60 Hz, 1300W	
Safety	4 fold safety system		
Compliance	CE marked, EU pressure directive 27/93/EG, Low voltage directive 2006/95/EG, Laboratory safety directive EN 61010-1 and -4, Electro magnetic compatibilty EN 61326, Workplace safety directive 2004/40/EG and 2008/46/EG		

# TECHNICAL DETAILS

	MULTICONTROL 2 TOUCH	CONNECT	CULTURE MEDIA
Key benefit	Freely programmable autoclave with documentation in the cloud	Freely programmable autoclave, with electrical lock and documentation in the cloud	Culture Media preparator with internal stirrer and process ports
Core use case	Sterilisation und materials testing in environments with high documentation requirements	Fully automatic autoclave for liquids sterilisation and documentation in the cloud	Preparation of up to 10 litres of culture media
Temperature control	Electronic		
Interface	HD touch screen controller		
Pressure control	Via temperature	Electronic pressure sensor	Electronic pressure sensor
Steam release	Automatic Steam release cock	Electrical valve	Electrical valve
Other key components	WiFi, Printer port, PC Port (USB)	WiFi, Bluetooth	Stirrer, 10 l media bucket, media sensor, process ports, WiFi, Bluetooth
Optional equipment	wire basket, printer, holder for for tablet controller	wire basket, media sensor, holder for tablet controller	wire basket, holder for tablet controller, port for additives, extraction port
Temperature range	Room temp. to 140° C, water bath mode with open lid		
Max working pressure	2,7 bar		
Chamber material	Stainless steel 304L		
Volume	18		
Diameter	18		
Rated voltage, Output	230 V AC, 50-60 Hz, 1300 W		
Safety	4 fold safety system		
Compliance	CE marked, EU pressure directive 27 / 93 / EG, Low voltage directive 2006 / 95 / EG, Laboratory safety directive EN 61010-1 and -4, Electro magnetic compatilbity EN 61326, Workplace safety directive 2004 / 40 / EG and 2008 / 46 / EG		

## CONTACT

---



### **CERTOCLAV**

Certoclav Sterilizer GmbH  
Georg Grinninger-Str. 37  
A-4050 Traun  
Phone: +43 (0)7229 689  
[www.certoclav.com](http://www.certoclav.com)

**Orders, technical customer service, custom-manufacture, dealer verification:**

Christina Bokesch, Phone: +43 (0)7229 689 26, Fax: +43 (0)7229 689 20  
Email: [bokesch@certoclav.com](mailto:bokesch@certoclav.com)

**Product Manager:** Klaus David, Phone: +43 (0)7229 689 28, Email: [david@certoclav.com](mailto:david@certoclav.com)

**Operating instructions, technical information and operation videos:** [www.certoclav.com](http://www.certoclav.com).